




# CONTAINMENT & SAMPLING T E C H N O L O G Y

- Multi-functional Glovebox system is integrated with the ProSys Recirculation Top Reactor Sampling System (RTRS)
- Multi-function glovebox that can be mounted either for liquids/solids charging directly into the sight-glass assembly of the RTRS or for the contained removal of samples from the RTRS
- Transfer in/out of the glovebox via Rapid Transfer Port
- Working chamber is nitrogen purged
- Internal safe change HEPA filters (removal via Rapid Transfer Port)
- Complete system can also be trolley mounted, when space on top of the reactor is limited

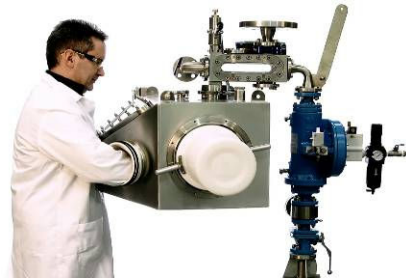

ProSys Sampling & Slurry Charging Gloveboxes

• info@prosys.ie •  www.prosys.ie •  +353 21 4613890



www.prosys.ie